

Patient Instructions


Pediatric Urine Collection

Your health care provider has requested a urine specimen.


In order to produce a meaningful laboratory test result, follow these instructions exactly.

Pediatric Urine Collector

Directions for Use:

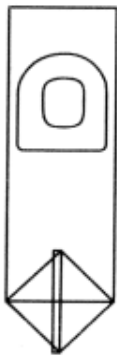

1. Separate the infant's legs. Make sure perineal and pubic areas are clean and dry.
2. Peel off backing paper, exposing adhesive.
3. For Female Use: Start bottom adhesion at skin flap separating rectum from vagina and work upward. Press firmly to skin surrounding the vagina with excess of bag toward feet.
4. For Male Use: Position penis naturally through bag operating with excess of bag toward the feet. Press adhesive firmly to entire perimeter without puckering.
5. Reapply diaper, if required, directly over bag.

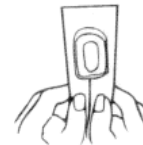

6. Remove and fold bag along long axis of the opening, thus folding adhesive surface in half. Press adhesive surfaces together firmly to assure a tight seal.

Labelling

1. Label the container with the child's name, date of birth, and the date of collection.

Packaging

1. Place the container in the bag provided and securely close the bag.
2. Fold any paperwork and place in the external pocket of the specimen transport bag. This ensures the specimen will not leak onto the paperwork.

Storage and Transport

1. Refrigerate and bring specimen to the laboratory **within one hour of collection**.